



- Program consists of an exe file and directory to hold dicom files





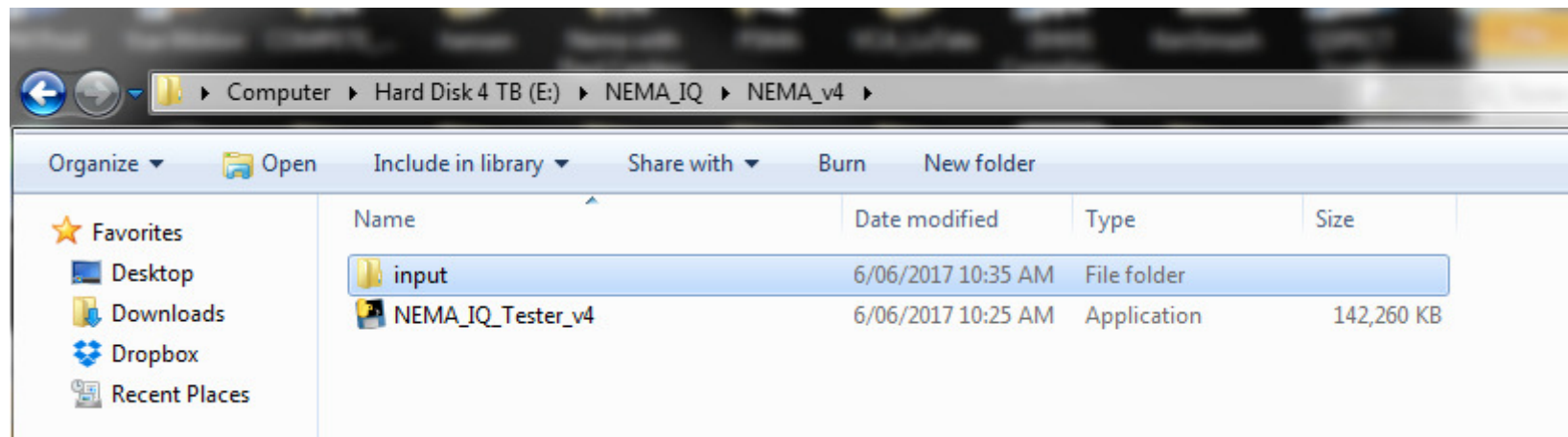

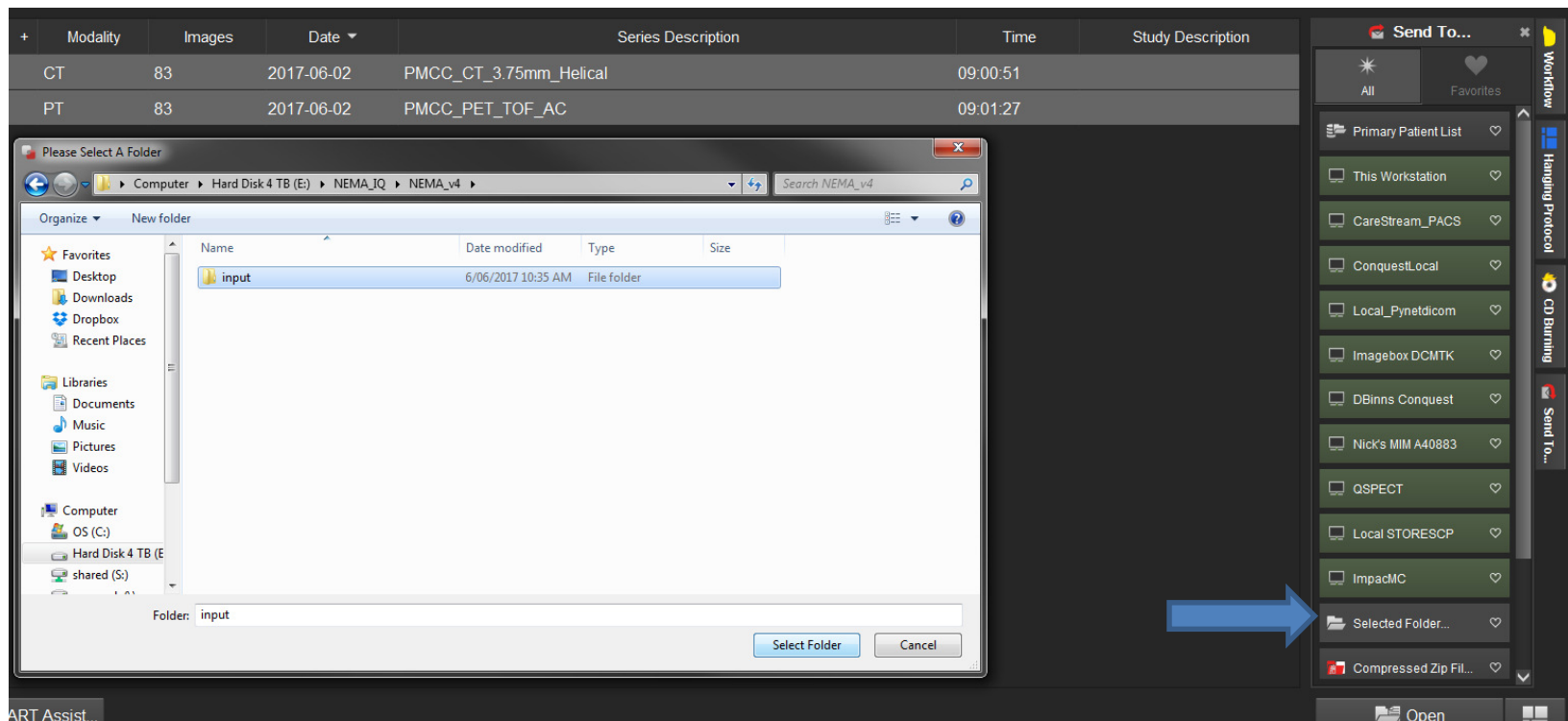

- Delete the contents of the input directory and copy your NEMA scan into the folder.
- Don't worry if CT is included or there are nested subdirectories. It will sort through and find the right scan



- Double-click the IQ tester program
- A bunch of text will pop up once it starts running

A screenshot of a Windows application window titled "E:\NEMA\_IQ\NEMA\_v4\NEMA\_IQ\_Tester\_v4.exe". The window displays a large grid of numbers arranged in rows and columns. Each row contains 20 numbers, starting from 23 and ending at 831. The numbers are organized into groups of four columns per row, with some rows having additional numbers or labels like "[1.", "2.", "3.", "4.", "5.", "6.", "7.", "8.", "9.", "10.", "11.", "12.", "13.", "14.", "15.", "16.", "17.", "18.", "19.", "20.", "21.", "22.". The grid is scrollable, as indicated by the vertical scrollbar on the right side. The background of the window is black, and the numbers are white.

- It should create a PDF report and some CSVs with the name of the scanner and the time of acquisition

 input	6/06/2017 10:39 AM	File folder	
 IQ_Data_D710_1_20170602 090127	6/06/2017 10:41 AM	Microsoft Excel C...	13 KB
 IQ_Results_D710_1_20170602 090127	6/06/2017 10:41 AM	Microsoft Excel C...	2 KB
 IQ_segmented_8_1_D710_1_20170602 090...	6/06/2017 10:41 AM	Adobe Acrobat D...	254 KB
 NEMA_IQ_Tester_v4	6/06/2017 10:25 AM	Application	142,260 KB
 qt.conf	6/06/2017 10:41 AM	CONF File	1 KB

- Background SUV is also included

